



natural-fill

FOR NATURAL BEAUTY

natural-meso

FOR NATURAL BEAUTY

STERILE SINGLE USE FILLER CANNULA / NEEDLE


natural-fill ?

Natural fill can be performed through a single access point, and the extended length of a cannula allows for the selection of the safest incision area. Also it can be bent around tissue structures and push away blood vessels and nerve, where a needle would simply pass through. Thus, cannula can reduce pain and shorten recovery time, making the procedure more comfortable for the patient. Furthermore, blunt cannula helps to reduce the risks of intravenous arterial embolism, enhancing patient safety.

- Less downtime
- Filling the solution
- Liposuction
- Safety, do not injure the vein
- Filling the several places with one cannula

Ultra Thin Wall Cannula (UTW)

In the injection process, some doctors have problems with injection due to strong viscosity, cohesiveness and gel hardness of solution. Therefore, we provide the Ultra Thin Wall Cannula (UTW) that have wider inner diameter for smooth injection and is more flexible than normal cannula.

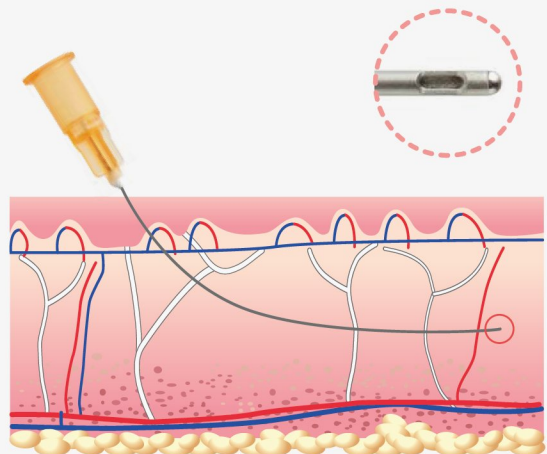

Anatomical Region	Cannula Recommendation
Periorbital	25, 27, 30G : 25 ~ 40mm
Forehead	25, 27G : 38 ~ 50mm
Brow	25, 27G : 25 ~ 50mm
Temporal, Middle third of the face, Lower third of the face	25, 27G : 38 ~ 50mm
Lips	25 - 30G : 25 ~ 40mm
Chin, Neck, Decollete, Hands	25 - 27G : 38 ~ 50mm

Gauge	Length
17G	70 / 90mm
18G	50 / 60 / 70 / 90mm
19G	50 / 60 / 70 / 90mm
20G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm
21G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm
22G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm
23G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm
24G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm
25G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm
26G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm
27G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm
29G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm
30G	25 / 30 / 35 / 40 / 50 / 60 / 70 / 90mm

natural-fill
FOR NATURAL BEAUTY

alfa

Flexible Cannula

Filler cannula is designed with a rounded tip and a flexible cannula which outstandingly smaller the possibility of injury to blood vessels contracted with other types of needles.



CE 1370

Graduation shows the injection volume and assists users' delicate control for the best procedure results.

Advanced cannula

Fine processing applied cannula provides outstanding uniformity, superb flexibility, and excellent performance.


natural-fill
FOR NATURAL BEAUTY
beta

Advanced Cannula with Graduation



Polycarbonate Hub

Polycarbonate is a thermoplastic which has a light yet hard property. It is classified as an engineering plastic with the features of excellent strength, heat resistance, impact resistance, and transparency.

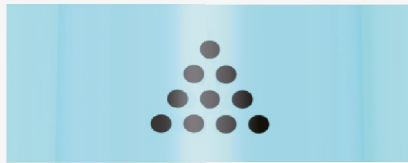

Indicator for the accurate control

Indicator is designed for ease to find the hole at the tip. This visually supports practitioners for the accurate control to give the superb treatment.



Graduation for the delicate injection

Graduation shows the injection volume and assists users' delicate control for the best procedure results.



UV adhesive

UV (ultraviolet rays) adhesive has advantages such as High adhesive strength and precision fixation, excellent chemical resistance, moisture resistance, and heat resistance, and its transparent characteristic makes it possible to realize an excellent design.

natural-fill
FOR NATURAL BEAUTY
combo

Advanced Cannula with Graduation + Guide Needle

*** Guide Needle : Puncture the skin in order to get into the intended space before using cannula.**

For smooth procedure, **Guide needle and natural-fill beta** are packed together. Whatever size of cannula you use, it is easy to puncture the skin by using the Guide needle packed in natural-fill combo.

Cannula		Guide Needle	
Gauge	Length	Gauge	Length
30G	25mm	27G	13mm
27G	40mm	26G	13mm
25G	50mm	23G	25mm
25G	40mm	23G	25mm
23G	50mm	21G	25mm
22G	70mm	21G	25mm
22G	50mm	21G	25mm
18G	50mm	-	-



CE 2292



natural-meso

FOR NATURAL BEAUTY



The thin wall needle enables less extrusion force and results in smoother injection. Meso therapy needle is effectively used for various ways such as injecting hyaluronic acid, local anesthetics, mesotherapy and placenta as well as high viscid solution.

Gauge	Length
27G	4 / 13mm
30G	4 / 6 / 8 / 10 / 13mm
31G	4 / 6 / 8 / 10 / 13mm
32G	4 / 6 / 8 / 10 / 13mm
33G	4 / 6 / 8 / 13mm
34G	4 / 6 / 8 / 13mm